

- across the body for approximately 60 seconds per region
- Only use the ELEEEELS percussive massage device on the body's soft tissue as desired. Do not use on the head or any hard and bone areas of the body. If any pain occurs during use, stop immediately.
  - Do not place any objects into the ventilation holes of the ELEEEELS percussive massage device.
  - Do not allow the device to come into contact with water.
  - Do not drop the device.
  - Only recharge with the supplied AC adapter.
  - Do not tamper or alter the device in any way.
  - Never leave the device operating or charging unattended



#### DISCLAIMER

PLEASE DO NOT USE THE ELEEEELS PERCUSSIVE MASSAGE DEVICE, OR ANY OTHER VIBRATION DEVICE WITHOUT FIRST OBTAINING APPROVAL FROM YOUR DOCTOR, SPECIALIST OR PHYSIOTHERAPIST, especially if you are pregnant, have diabetes with complications (such as neuropathy or retinal damage), wear a pacemaker, have had recent surgery, suffer from epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, have had recent joint replacements or IUD'S, metal pins or plates, or you have any general concerns about your physical health. Frail Individuals and children should be accompanied by an adult when using any percussive massage devices.

#### WARRANTY

IF WITHIN ONE YEAR FROM THE DATE OF PURCHASE, THIS PRODUCT FAILS DUE TO A DEFECT IN MATERIAL OR WORKMANSHIP, ELEEEELS WILL REPAIR OR REPLACE THE PRODUCT, OR NECESSARY COMPONENTS, FREE OF CHARGE

#### THIS WARRANTY EXCLUDES:

- (a) damage caused by accident, abuse, mishandling, or transport;
- (b) units subjected to unauthorized repair;
- (c) units not used in accordance with ELEEEELS care instructions;
- (d) damage exceeding the cost of the product;
- (e) deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) failure to show the dated proof of purchase



ELEEEELS

Unit 5, 15/F, Global Gateway Tower,  
63 Wing Hong Street, Cheung Sha Wan,  
Kowloon, Hong Kong  
[www.eleeeels.com](http://www.eleeeels.com)

X1TE2



# X1T

Percussive  
Massage  
Device

## Instruction Manual




ELEEEELS

We have developed the ELEELLS percussive massage device to relieve pain and promote circulation, which helps to accelerate recovery and prevents muscle damage.

**CORDLESS/RECHARGEABLE**  
24 V Lithium-ion Battery

**LIGHTWEIGHT DESIGN**  
2.5lbs / 1.1kg

## SPECIFICATIONS

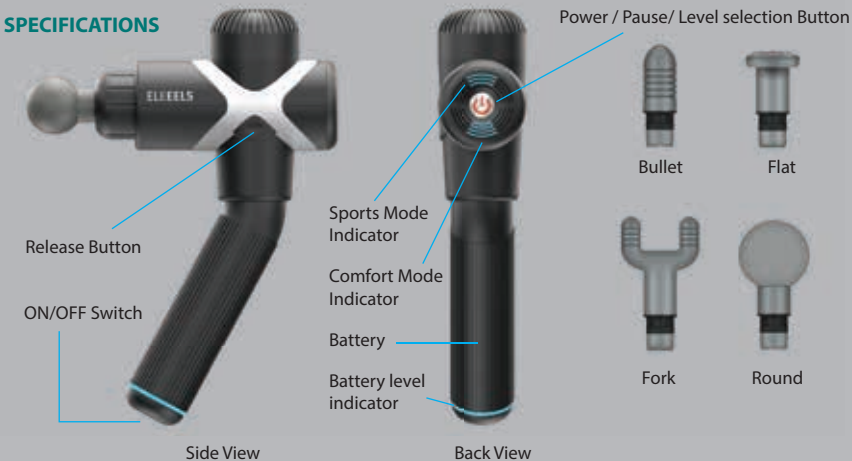

- 1 Powerful Brushless High-Torque Motor
- 2 Removable Battery Release
- 3 Rechargeable Lithium-Ion Battery (3+ Hrs of Use Per Charge)
- 4 6 Speed Settings for Optimal Pressure
- 5 Cordless and Lightweight (2.5 lbs)
- 6 ON / OFF Power Switch
- 7 Interchangeable Applicators (x4)
- 8 ELEELLS universal adaptor included

Vibration Frequency:  
Comfort Mode: 900rpm, 1200rpm, 1500rpm  
Sports Mode: 1800rpm, 2400rpm, 3200rpm

Weight: 2.5 lbs / 1.1 kg  
Charge: 110-240V 50/60Hz 28V, DC-1A

## USAGE INSTRUCTIONS

PLEASE READ PRIOR TO USE

### CHARGING

Fully charge the battery for up to 6 hours before first use, either attached or detached from the device

To charge the battery:

- 1) Turn off the main power switch located at the bottom of the handle and connect to the supplied adaptor.
- 2) A green light will flash at the bottom of the handle when charging has begun
- 3) When the green light has stopped flashing, the device has been fully charged

### OPERATING

- Attach the battery/handle into the base of the device unit.
- Press the main power button at the bottom of the handle to "ON".
- Press the power button located at the rear of the device.
- The massage gun will stop operation after a 10 minute session, simply press the power button to start a new one.

Sports Mode level 1-3:

Press the bottom power button once to turn device on and into level 1 of sports mode. Pressing the main power button can adjust the speed to level 2 and 3 respectively. The blue LED lights will indicate the corresponding power level.

Comfort Mode level 1-3

Hold the main power button for 2 seconds and release to turn device on and into level 1 of comfort mode. Further pressing the main power button can adjust the speed to level 2 and 3 respectively. The blue LED lights will indicate the corresponding power level.

- Massage the desired body part at the desired speed by applying pressure as needed

- To remove the battery, press the battery release button with your thumb. Firmly pull downwards on the battery with the other hand
- The circular light indicator on the battery:
  - Red: 0-25%
  - Blue: 25-75%
  - Green: 75- 100%

### HOW DO I ADJUST THE HANDLE'S ANGLE?

Twist the battery (handle) 180 degrees. A 'click' can be heard when this is completed

### MAINTENANCE, CLEANING & STORAGE

- Use a slightly damp towel to clean the device and a soft cloth to dry it.
- Do not clean the device with chemicals or alcohol.
- Only use the ELEELLS supplied AC adaptor.

### WARNINGS & CAUTION

- Do not immerse in water. Keep away from liquids or heat sources. Keep ventilation ports free from dust and debris.
- Do not remove screws or attempt to disassemble.
- Do not operate continuously for more than one hour. Allow device to rest for 10 minutes before re-using.
- For adult use only. Do not use if injured. Consult your doctor before using this product.

### SAFETY INSTRUCTIONS

DANGER: TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

- For adult use only
- Use only on a dry clean surface of the body over clothing by lightly pressing and moving