X4

Percussive
Massage
Device


Instruction Manual

CORDLESS / RECHARGEABLE
24V Lithium-ion Battery

LIGHTWEIGHT DESIGN
2.4 lbs. / 1.1 kg


OPERATING

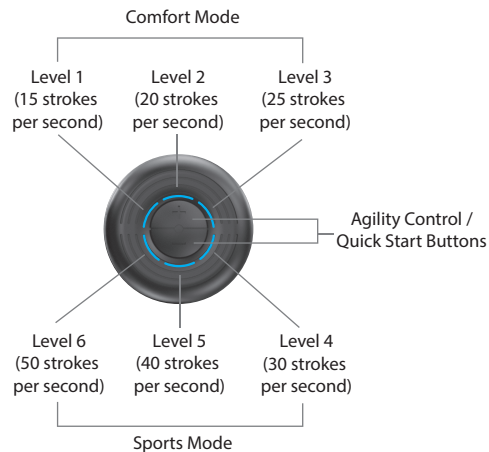- Switch on the main power button at the bottom of the handle.
- Sports Mode Quick Start: Press the "+" agility control button once for level 4.
- Comfort Mode Quick Start: Press the "-" agility control button once for level 3.
- Press the "+" or "-" agility control button repeatedly to adjust the speed or pause the device during operation. The blue LED lights indicate the corresponding speed level.
- Massage the desired body part at your preferred speed by applying pressure as needed.
- When the massage device stops operating after a 10-minute session, press the agility control button to start a new session.
- Switch off the main power button after operation.

We developed the ELEEEELS percussive massage device to accelerate recovery and prevent muscle damage in order to relieve pain and promote circulation.

ELEEEELS X4

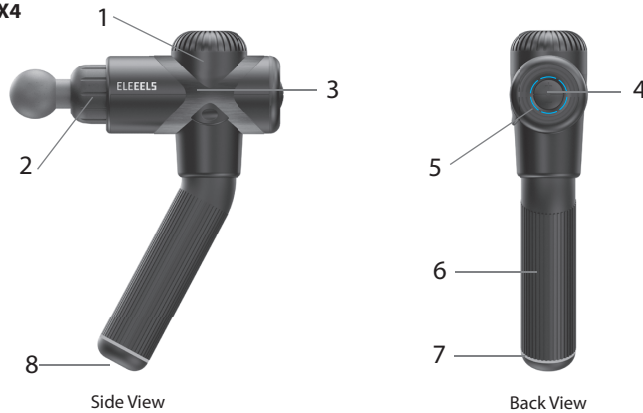- 1) High-torque, high-speed brushless motor (delivers up to 50 strokes per second)
- 2) Xtreme Power™ Engineering (13 mm amplitude deep muscle therapy)
- 3) Eterno Quiet™ Motion Technology
- 4) +/- 6 speed agility control for optimal pressure
- 5) Speed level indicator
- 6) Rechargeable lithium-on battery (up to 10 hours of use per charge)
- 7) Battery level indicator
- 8) ON/OFF switch

Vibration Frequency:
Comfort Mode: Level 1 (15 strokes per second), Level 2 (20 strokes per second), Level 3 (25 strokes per second)
Sports Mode: Level 4 (30 strokes per second), Level 5 (40 strokes per second), Level 6 (50 strokes per second)

Weight: 2.4 lbs. / 1.1 kg
Charge: 110-240V 50/60Hz 28V, DC-1A

HOW DO I ADJUST THE ANGLE OF THE HANDLE?

Twist the battery (handle) 180 degrees until you hear a 'click'.

MAINTENANCE & CLEANING

- Use a slightly damp towel to clean the device and a soft cloth to dry it.
- Do not clean the device with chemicals or alcohol.
- Only use the AC adapter supplied by ELEEEELS for recharging.

WARNINGS & CAUTIONS


- Do not immerse in water. Keep away from liquids or heat sources. Keep ventilation ports free from dust and debris.
- Do not remove screws or attempt to disassemble.
- Do not operate continuously for more than one hour. Allow device to rest for 10 minutes before re-using.
- For adult use only. Do not use if injured. Consult your doctor before using this product.

SAFETY INSTRUCTIONS

DANGER: TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

- Apply to a clean, dry clothed area of the body by lightly pressing and moving across the body for approximately 60 seconds per area.
- Only use the device on the body's soft tissue as desired. Do not use on head or any hard, bony area of the body. Discontinue use immediately if experiencing any pain.
- Do not drop the device.
- Do not tamper or alter the device in any way.
- Never leave the device unattended during operating or charging.

5 Interchangeable Applicators

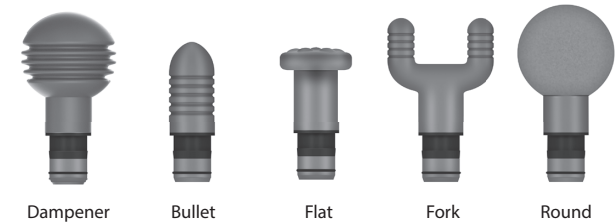USAGE INSTRUCTIONS

PLEASE READ PRIOR TO USE

CHARGING

Charge the battery for up to 6 hours before first use to make sure it is fully charged. To charge the battery:

- 1) Switch off main power button at the bottom of the handle, and connect the device to the supplied adapter.
- 2) A green light will flash at the bottom of the handle when the charging begins.
- 3) When the green light stops flashing, the device is fully charged.
- 4) The circular light power indicator shows the battery level.
 - a) Red: 0-25%
 - b) Blue: 25-75%
 - c) Green: 75- 100%

DISCLAIMER

Please do not use the ELEEEELS percussive massage device, or any massage device without first obtaining approval from your doctor if you have any of the following conditions: pregnancy, neuropathy, retinal damage, pacemaker user, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements, surgical pins or plates, or you have any concerns about your physical health.

WARRANTY

IF WITHIN ONE YEAR FROM THE DATE OF PURCHASE, THIS PRODUCT FAILS DUE TO A DEFECT IN THE MATERIAL OR WORKMANSHIP, ELEEEELS WILL REPAIR OR REPLACE THE PRODUCT; OR FAILED COMPONENT(S), FREE OF CHARGE

THIS WARRANTY EXCLUDES:

- (a) Damage caused by accident, abuse, or mishandling during transport;
- (b) Damage arising from use of units from unauthorized repair;
- (c) Malfunction of units that are not used in accordance with ELEEEELS care instructions;
- (d) Damage for which the cost of repair exceeds the cost of the product;
- (e) Deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) Any product for which a dated proof of purchase cannot be provided.

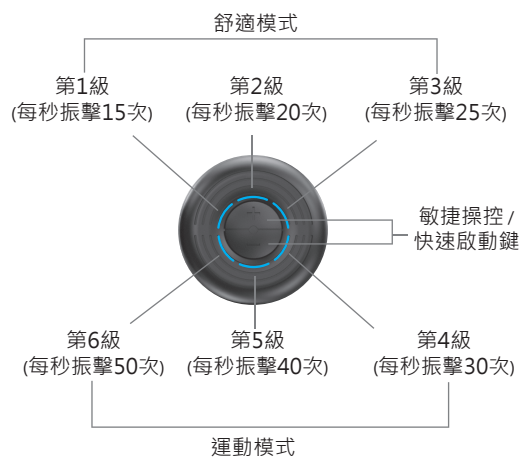
Please visit below link for instruction manual in other languages:
<http://www.eleeeels.com/im>

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